

## Sienna COVID-19 Antigen Rapid Test Kit

## Results in 10 minutes using a Nasopharyngeal swab

#### **Test Details**

The Sienna COVID-19 Antigen Rapid Test Cassette is a rapid chromatographic immunoassay for the qualitative detection of COVID-19 antigen in Nasopharyngeal swab. The identification is based on the monoclonal antibodies specific for the Nucleocapsid protein of SARS-CoV-2. It is intended for professional in-vitro diagnostic use only.

#### Accurate Results

Sensitivity: 95.0%

Specificity: 99.2%

Accuracy: 98.8%

#### **Unmatched Convenience**

Visual Results: Easy to interpret.

Fast Results: 10-minute test time.

Individual Buffer Vials and Sterile Swabs:

Helps in testing multiple individuals simultaneously.






### Sienna COVID-19 Antigen Rapid Test Kit

During the testing process, the extracted specimen reacts with the antibody of the SARS-CoV-2 nucleocapsid protein that are coated onto particles. The mixture migrates into the membrane to react with the antibody for the SARS-CoV-2 nucleocapsid protein on the membrane and will generate one colored line in the test regions. The presence of this colored line in the test regions indicates a positive result. A colored line will always appear in the control region to serve as a procedural control if the test has been performed correctly.

#### Test Procedure

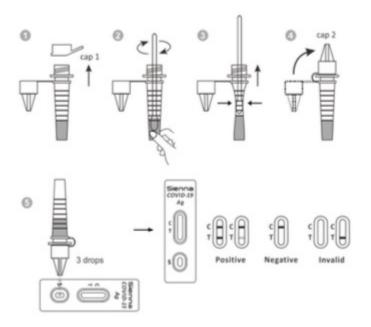

#### Kit Contains:

- (1) IFU
- (25) Individually pouched Cassettes
- (1) Workstation
- (25) Buffer Tubes
- (25) Sterile Swabs



# **Order and Shipping Information**

Cat. No. Product: QTY Registration Status

102242 Sienna COVID-19 Antigen Rapid Test 25 Tests CE

Contents: EUA Authorized May 2021

(1) IFU
(25) Individually pouched Cassettes
(1) Workstation
(25) Buffer Tubes
(25) Sterile Swabs

Item	Quantity 25 Tests	Length		Width		Height		Box Weight	
		5	in	9.35	in	2.66	in	0.57	Ibs
Box		12.25	cm	23.75	cm	6.75	cm	0.26	kg
Carton	49 Boxes	24	in	16	in	20	in	46	lbs
	(1,225 test)	60	cm	40	cm	50	cm	20.8	kg
Pallet	8 Cartons	48	in	42	in	72	in	370	lbs
- " -	0.000	100		105	cm.	180	cm	170	kg
together = 20	9,800 eight of Produc 0 kg Air Passer	ger Pla	ne Car	go Only.		er packagi	ng (25kg)		
NOTE: Total w together = 20 Pallet	eight of Produc 0 kg Air Passer 12 Cartons	t Appro ger Plan 48	x. (200k ne Carg in	go Only.	in	72	ng (25kg) in	550	lbs
NOTE: Total w together = 20 Pallet Option 2 NOTE: Total w	eight of Produc 0 kg Air Passer	48 120 1265kg	in cm	go Only. 42 105	in cm	er packagi	ng (25kg)		lbs
NOTE: Total w together = 20 Pallet Option 2 NOTE: Total w	eight of Produc 0 kg Air Passer 12 Cartons 14,700 eight of Produc	48 120 1265kg	in cm	go Only. 42 105	in cm	72	ng (25kg) in	550	lbs kg



# Sienna Rapid Antigen Test Purchase Process SOP



1: Purchase order made out to Adaptiv Biomed, Authorized US Distributor.

#### Also include:

- Order Number that we will assign.
- Any Sales Tax Exemption Form or Resale Certificate
- CLIA Number or CLIA Waiver (assistance can be provided in obtaining CLIA Waiver.
- NPI if doctor or medical practice
- 2: We will email confirmation of receiving PO, provide a Sales Order and propose delivery schedule.
- 3: Customer wires funds order to BLVD Medical. If financing, financing or terms are finalized.
- 4: Drop ships direct from US facility to Customer's desired delivery destination. International orders ship from Europe to desired destination.
  - Shipment confirmation and tracking provided.

